

MoerLab

T-Connect™ Tele-Coil Neckloop

Use Auracast™ audio
with T-Coil hearing aids






- T-Coil compatibility
- Magnetic induction
- 3.5mm Plug-and-play



System Connection Diagram

Auracast Transmitter ➔ MoerDuo Plus™ Receiver + T-Connect™ ➔ T-Coil Hearing Aids

- MoerDuo Plus™ receives the Auracast™ audio stream and delivers it to the T-Connect.™
- The T-Connect™ uses magnetic induction to relay the audio to T-Coil hearing aids.

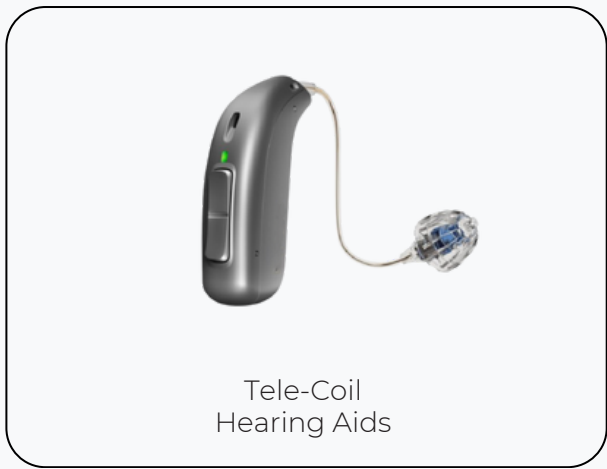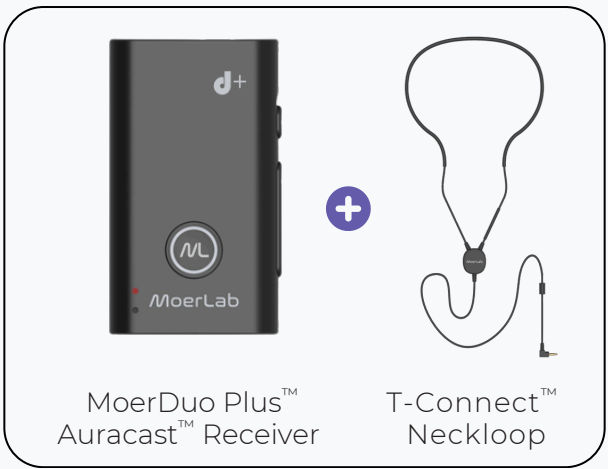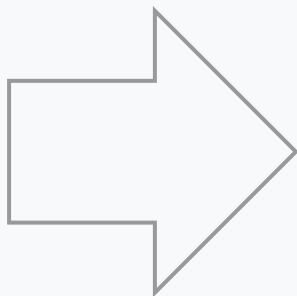


Key Features

- Moderate cable length
- Built-in safety release
- Low-impedance output
- Works with T-Coil hearing aids




Comfortable neck strap with built-in safety release

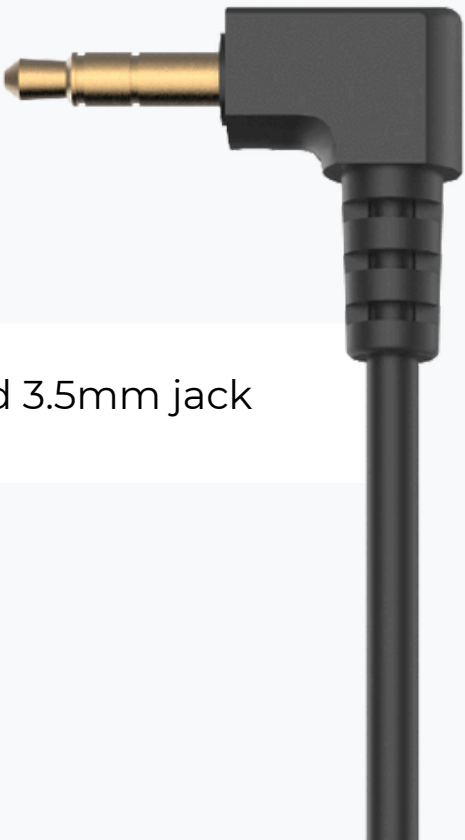

Angled 3.5mm jack



T-Connect™ connects to MoerDuo Plus™ via its angled 3.5mm audio plug

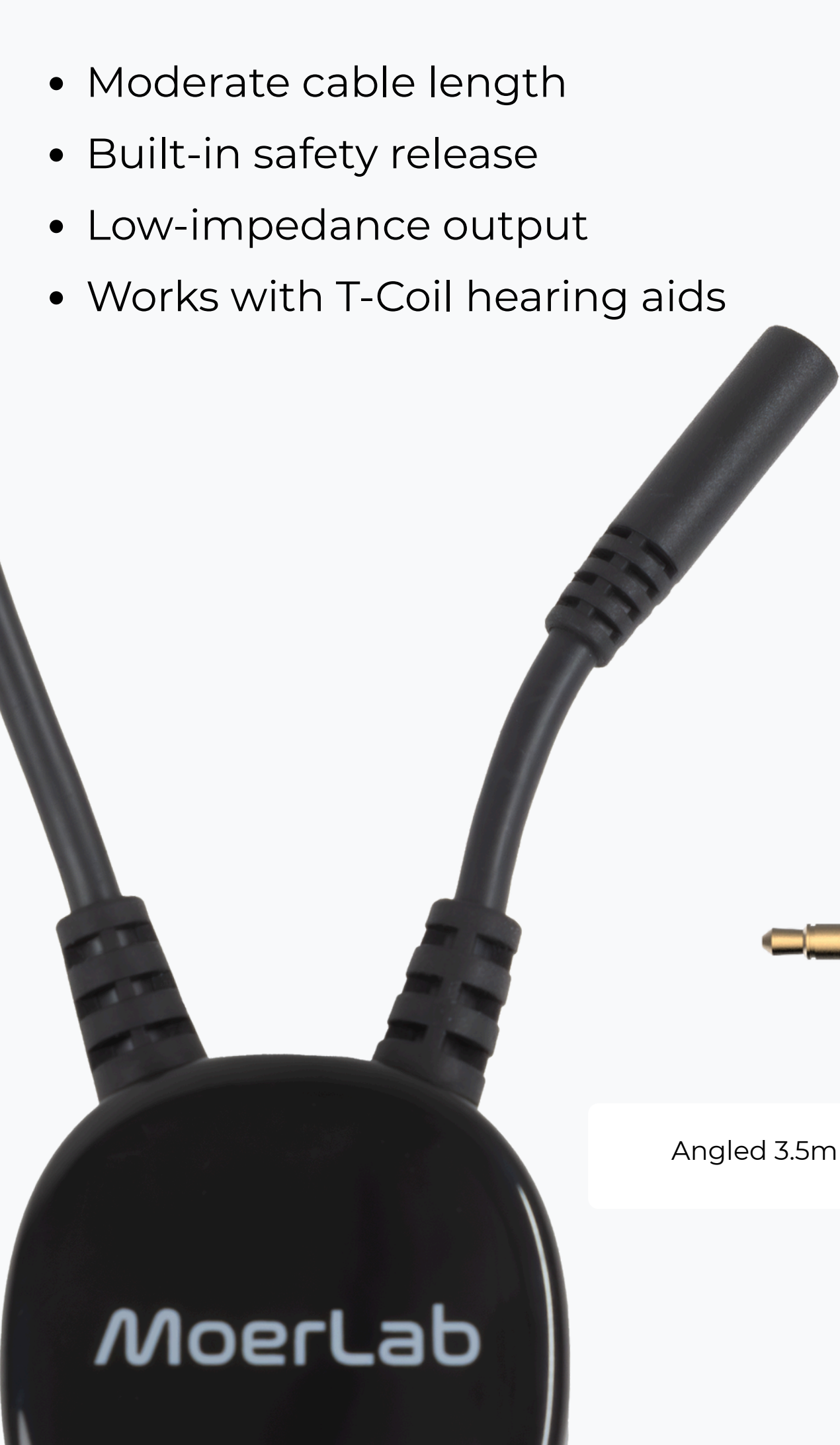

Experience clear audio via your T-Connect™ Tele-Coil Neckloop



Specifications

Product Name	T-Connect™ Tele-Coil Neckloop
Model	ATC01
Weight	70 g
Cable Length	121 cm
Connector	Angled 3.5 mm stereo plug (gold-plated)
Impedance	2Ω

